

# MagPurix<sup>®</sup> Blood DNA Extraction Kit 200 (ZP02001)

## Instructions for Use (Handbook)



Version: 2.1



48



**For *in vitro* diagnostic use**



**ZINEXTS LIFE SCIENCE CORP.**

16F., No. 93, Sec. 1, Xintai 5th Rd.,  
Xizhi Dist., New Taipei City 221416,  
Taiwan (R.O.C.)

**Obelis s.a.**

Bd Général Wahis 53  
1030 Brussels Belgium  
Tel: +(32) 2 732-59-54  
Fax: +(32) 2 732-60-03  
mail@obelis.net



Read and follow these Instructions for Use prior to using this product. The latest revision of this document can be found at [www.zinexts.com](http://www.zinexts.com)

# Contents

<b>Intended Use</b> .....	<b>3</b>
<b>Introduction</b> .....	<b>3</b>
<b>Kit Contents and Storage</b> .....	<b>4</b>
<b>Materials Required But Not Provided</b> .....	<b>5</b>
<b>Warnings and Precautions</b> .....	<b>5</b>
<b>Purification Principle</b> .....	<b>6</b>
<b>Before Starting</b> .....	<b>7</b>
<b>Preparation of sample materials</b> .....	<b>7</b>
<b>Procedure of MagPurix System</b> .....	<b>8</b>
<b>Purification Protocol - MagPurix<sup>®</sup> series</b> .....	<b>9</b>
<b>Purification Protocol - MagPurix<sup>®</sup> EVO series</b> .....	<b>9</b>
<b>Troubleshooting</b> .....	<b>11</b>
<b>Related Products</b> .....	<b>12</b>
<b>References</b> .....	<b>12</b>
<b>Symbols</b> .....	<b>13</b>
<b>Limited Product Warranty</b> .....	<b>14</b>
<b>Revision History</b> .....	<b>14</b>

## Intended Use

The MagPurix<sup>®</sup> Blood DNA Extraction Kit 200 provides a complete set of reagents and consumables for the automated purification of genomic DNA (total nucleic acids) from mammalian whole blood, buffy coat and leukocyte concentration with MagPurix system.

The product is intended to be used by professional users, such as technicians and physicians who are trained in molecular biology techniques.

## Introduction

Product Name	MagPurix <sup>®</sup> Blood DNA Extraction Kit 200
Catalogue Number	ZP02001
Product Overview	The MagPurix <sup>®</sup> Blood DNA Extraction Kit 200 is designed to extract genomic DNA from mammalian whole blood, buffy coat and leukocyte concentration. The kit uses unique magnetic ZiBeads <sup>®</sup> technology and in combination with MagPurix <sup>®</sup> series automatic instruments, achieved superior product quality, consistent and high product yield and reproducible results. The final product is suitable for a wide range of diagnostic and research applications, such as sequencing, genotyping, qPCR, ddPCR and NGS assays.
Applicable Instrument Model	All MagPurix <sup>®</sup> Instrument
Display Protocol Name on The Instrument	2001 BLOOD DNA 2001 BLOOD DNA RAPID (EVO only)
Applicable Instrument Firmware	Check and download the latest firmware from <a href="http://www.zinexts.com">www.zinexts.com</a>
Processing Time	MagPurix <sup>®</sup> 12 series 55-75 minutes MagPurix <sup>®</sup> 24 series 60-80 minutes MagPurix <sup>®</sup> EVO series 50-65 minutes ( <b>RAPID</b> : 30 minutes)


## Kit Contents and Storage

Shipping and Storage	The kit is shipped at room temperature. Upon receipt, store the kit at room temperature. All kit components are stable when stored properly until the expiration date shown on the kit box.	
Kit Content	The components supplied in the kit are listed below. Sufficient reagents are supplied to perform 48 purifications.	
	Contents	Amount
	<b>1</b> Reagent Cartridge	48 pcs (6x8)
	<b>2</b> Reaction Chamber	48 pcs (6x8)
	<b>3</b> Tip Holder	48 pcs (6x8)
	<b>4</b> Piercing Pin	50 pcs
	<b>5</b> Filter tip	50 pcs
	<b>6</b> Sample Tube (2 ml)	50 pcs
	<b>7</b> Elution Tube (1.5 ml)	50 pcs
	Barcode sticker (EVO only)	50 pcs

Reagent  
Cartridge  
Contents

Each Reagent Cartridge has 10 positions with 10 sealed well.  
Positions 1-10 contain wells filled reagents for this protocol.

Reagent	Well No.
Proteinase K Solution	1
Lysis Buffer 2	2
Binding Buffer 8	3
Magnetic Bead Solution	4
Washing Buffer 1	5
Washing Buffer A	6
Washing Buffer B	7
Elution Buffer 1	8
Elution Buffer 2	9
Empty	10

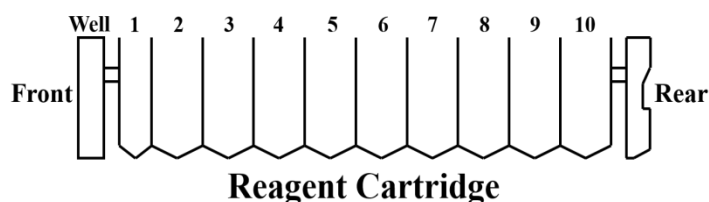

## Materials Required But Not Provided

The following general laboratory equipment and consumables are required to perform the extraction. All laboratory equipment should be installed, calibrated, operated, and maintained according to the manufacturer's recommendations. The following tables display required and special equipment along with the list of consumables.

Item
MagPurix <sup>®</sup> series instrument
1.5 or 2.0 ml microcentrifuge tubes
Pipettes and filter tips
Phosphate-buffered saline (PBS, may be required for diluting samples)
Optional: DNase-free RNase A (to minimize RNA content)

## Warnings and Precautions

For *in vitro* diagnostic use only. Read all the instructions carefully before using the kit. Use of this product should be limited to trained personnel in the techniques of DNA purification. Strict compliance with the user manual is required for optimal results. Attention should be paid to expiration dates printed on the box and labels of all components. Do not use a kit after its expiration date.


When working with chemicals, always wear a suitable lab coat, disposable gloves, and protective goggles. For more information, please consult the appropriate safety data sheets (SDSs). These are available online in convenient and compact PDF format at <http://www.zinexts.com/download.php?lang=en&tb=1&cid=7> where you can find, view, and print the SDSs for each kit and kit component.

Please report any serious incident occurred in relation to the device to your local representative/agent or the manufacturer, and to the competent authority of your country/state.



CAUTION: DO NOT add bleach or acidic solutions directly to the sample preparation waste.

# Purification Principle

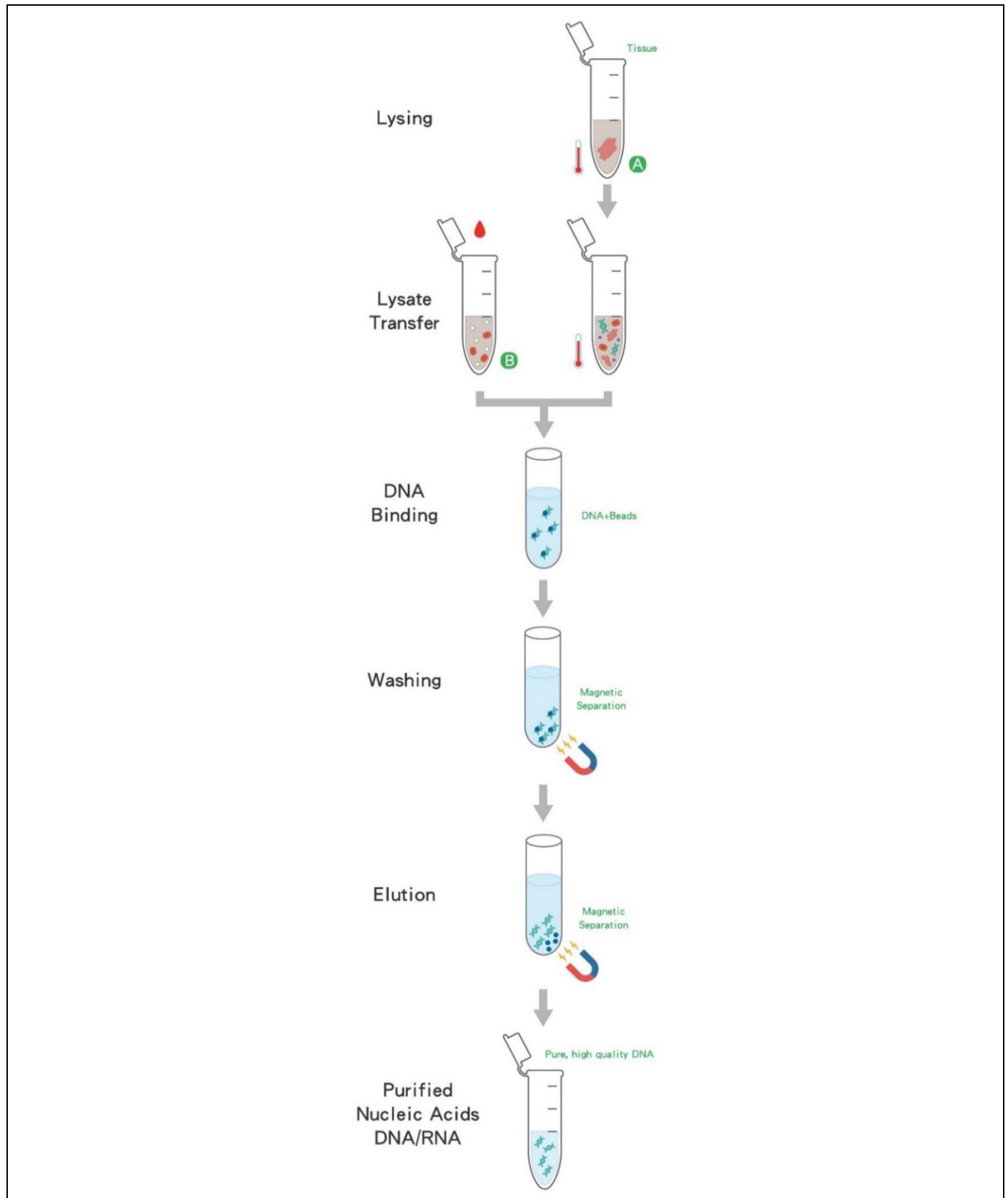

- A** Perform certain pretreatment process before extraction.
- B** Transfer sample to extraction directly.

# Before Starting

## Preparation of sample materials

The purification procedure is optimized for the use of 100-400 µl whole blood, buffy coat and leukocyte concentration.

Mammalian Whole Blood	<ol style="list-style-type: none"><li>Gently invert the blood collection tube.</li><li>Dispense 100-400 µl blood sample into each Sample Tube.</li><li>If the sample volume is lower than 100 µl, please complete the volume with appropriate amount of PBS.</li><li>If the white blood cell number of the sample is more than <math>2 \times 10^4</math> cells/µl, please dilute the sample with PBS into appropriate concentration.</li></ol>
Buffy Coat	<ol style="list-style-type: none"><li>Dilute the buffy coat sample to 20 times of its volume with 1X PBS and mix gently.</li><li>Adjust the sample volume to ensure that nucleated cells do not exceed the concentration of <math>2 \times 10^4</math> cells per µl.</li><li>Dispense 100-400 µl buffy coat sample into each Sample Tube.</li></ol>
Leukocyte Concentration	<ol style="list-style-type: none"><li>If sample contain less than <math>1 \times 10^3</math> leukocytes/µl, please concentrate the blood cells by centrifuge at <math>300 \times g</math> for 15 minutes at 4°C and use the concentrated leukocytes for extraction procedure.</li><li>Discard plasma</li><li>Dispense 100-400 µl leukocyte concentrate (no more than <math>5 \times 10^6</math> cells) into each Sample Tube.</li></ol>

### Note:

This kit is used to work with fresh or frozen blood samples collected in tubes which contain common anti-coagulants like EDTA, heparin\* and citrate. (\*EDTA is recommended to use as anticoagulation agent, while heparin has inhibitory effects on nucleic acid amplification reaction.)

Using fresh whole blood sample (within 1 week, stored at 4-8°C) for extraction is recommended. Total nucleic acid yield and quality will decrease with time or after multiple thawing. For longer storage time, whole blood should be frozen at -20°C or lower and avoid freeze-thaw cycles.

For use of the concentrated buffy coat (purified and free of blood cells), the MagPurix<sup>®</sup> Cultured Cell DNA Extraction Kit (ZP02005) is recommended.

This protocol was established for isolating DNA from whole blood of healthy individual; affected individual or drug-treated individual (e.g., Patients affected by leukemia or infection disease) may have abnormal blood quality and influence the nucleic acid extraction procedure. **If the whole blood sample is granulocyte-rich (white blood cell number more than  $2 \times 10^4$  cells/µl), dilute blood sample with PBS or extract the DNA by the MagPurix<sup>®</sup> Blood DNA Extraction Kit 1200 (ZP02002).**

The final eluate contains total nucleic acid (DNA and RNA). RNA is not the major product in this kit (about 10%) and would degrade soon. If the RNA-free product is needed, please add RNase to treat the eluate. (For RNase treatment, follow the manufacturer instructions of the kit used in your lab.)


The suggested starting material and elution volume range for each nucleic acid extraction.

Sample type	Starting material per sample	Elution Volume
Mammalian Whole Blood	100-400 µl whole blood NOTE: When the WBC number is more than $2 \times 10^4$ cells/µl, sample dilution with PBS into appropriate concentration is recommended.	100-300 µl
Buffy Coat*	100-400 µl diluted buffy coat* *Must dilute the buffy coat to at least 20 times of its volume with 1X PBS.	
Leukocyte Concentration**	100-400 µl leukocyte concentrate **The concentration contains no more than $5 \times 10^6$ cells. Samples with low leukocyte count (less than $1 \times 10^3$ cells/µl), please concentrate the blood cells by centrifuging at $300 \times g$ for 15 minutes at 4°C before DNA extraction.	

⚠ Treat all samples as potentially infectious, follow the local safety regulation

## Procedure of MagPurix System

### Workflow of MagPurix operation

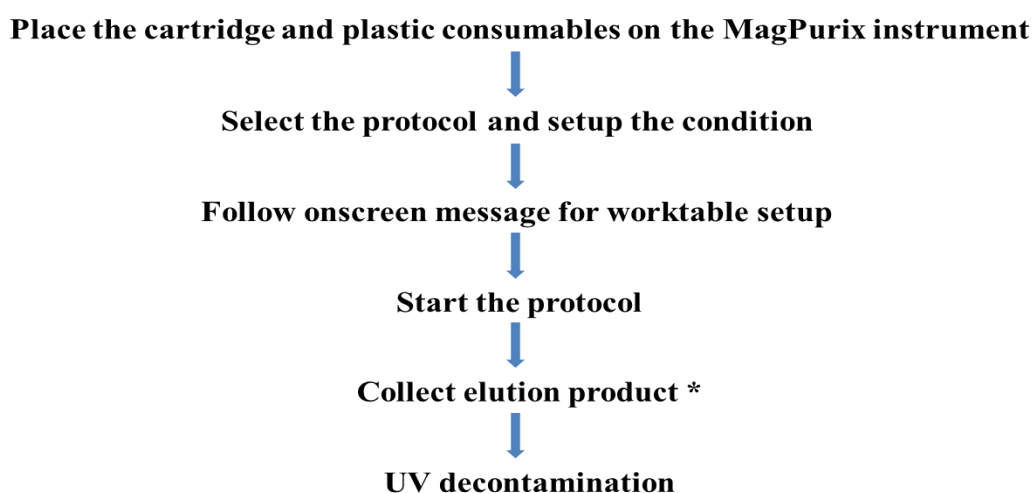

\* Output the bench record (option)

Note: Perform all steps at room temperature (20-25°C) unless otherwise notified.





## Purification Protocol - MagPurix<sup>®</sup> series

1	Turn on the Instrument	a. Turn on the power switch and wait for the screen to turn on.
2	Load new Consumable(s) and Cartridge(s)	a. Open the door and remove the Sample Rack from the instrument. b. Load <b>1</b> Reagent Cartridge, and all plastic disposables ( <b>2</b> Reaction Chamber, <b>3</b> Tip Holder, <b>4</b> Piercing Pins, <b>5</b> Filtered Tips and other components presented in the kit intended to use). c. Place <b>6</b> Sample Tubes and <b>7</b> Elution Tubes into the Sample Rack.
3	Load the Samples	a. Transfer appropriate volume of sample into each Sample Tube on the Sample Rack. b. Put the Sample Rack back into the instrument and close the door.
4	Program Set up	a. Scan the protocol barcodes to select the purification protocol, sample volume and elution volume.
5	Start Extraction	a. Follow the instructions displayed on the screen to double-check the operating steps being completed before program running. b. Press “ <b>ENTER</b> ” to start the experiment. Instrument will run the protocol program automatically until the whole process is completed. c. At the end of the run (approximately 12 series <b>55-75 minutes</b> , 24 series <b>60-80 minutes</b> ), instrument alarms briefly.
6	Collect the Elution tubes	a. Open the instrument door. b. Collect the elution tubes containing the purified nucleic acids. c. The purified nucleic acids are ready for immediate use. Store the purified nucleic acids at 4°C (short-term, less than 10 days) or aliquot and store at -70°C (long-term) before performing downstream analysis. d. Discard the used cartridges and all plastic consumables into biohazard waste. <b>*Do not reuse the cartridges.</b> e. If you are not using the instrument immediately, place the Sample Rack back to the workplace, close the instrument door and press “Start” button for 2 seconds to enter sleep mode. Moreover, if the instrument will not be used in a long time, please turn off the power switch.

## Purification Protocol - MagPurix<sup>®</sup> EVO series

1	Turn on the Instrument	a. Turn on the power switch and wait for the screen to turn on. b. Login the instrument and enter the Home Page.
2	Load new Consumable(s) and Cartridge(s)	a. Open the door and remove the Sample Rack from the instrument. b. Open the Tip-Holder Lid. c. Load <b>1</b> Reagent Cartridge and all plastic disposables ( <b>2</b> Reaction Chamber, <b>3</b> Tip Holder, <b>4</b> Piercing Pins, <b>5</b> Filtered Tips and other components presented in the kit intended to use). d. Close the Tip-Holder Lid.

- 
- e. Paste the barcode stickers on Elution Tubes.  
f. Place **6** Sample Tubes and **7** Elution Tubes into the Sample Rack.
- 
- 3** Load the Samples
- a. Transfer appropriate volume of sample into each Sample Tube on the Sample Rack.  
b. Put the Sample Rack back into the instrument and close the door.
- 
- 4** Program Set up
- a. Select the appropriate protocol program on the instrument. Press **NEXT**.  
b. Select the appropriate Sample Volume and Elution Volume and press **NEXT**.  
c. Press the number button to select the right Sample Numbers.  
d. Scan/Edit each primary Sample ID directly. After finished, press **NEXT**.  
e. Scan/Edit each Elution Tube ID directly. After finished, press **NEXT**.  
f. Scan Reagent Cartridge Barcode. Press **NEXT**.  
*\*If the cartridge is expired, the next step cannot be performed.*  
g. Follow the instructions on the screen to double-check the operating steps being completed before running the program. Press **NEXT**.
- 
- 5** Start Extraction
- a. Check "**PROGRAM CONFIRMATION**" on the screen.  
b. Press "**START**" to start the experiment. Instrument will run the protocol program automatically until the whole process is completed.  
c. At the end of the run (approximately **50-65 minutes**) (**RAPID: 30 minutes**), instrument alarms briefly and the screen indicates "**PROGRAM FINISH**".  
d. If you want to perform the same protocol, press "**RERUN**" to perform the same experiment. If you do not need to re-run the experiment, press the function button " **HOME**" to exist the experiment mode.
- 
- 6** Collect the Elution tubes
- a. Open the instrument door.  
b. Collect the elution tubes containing the purified nucleic acids.  
c. The purified nucleic acids are ready for immediate use. Store the purified nucleic acids at 4°C (short-term, less than 10 days) or aliquot and store at -70°C (long-term) before performing downstream analysis.  
d. Discard the used cartridges and all plastic consumables into biohazard waste. *\*Do not reuse the cartridges.*  
e. If you are not using the instrument immediately, please put the Sample Rack back into the instrument, close the instrument door, and press the " **POWER**" function button to enter sleep mode. If the instrument will not be used in a long time, please turn off the power switch.
-

# Troubleshooting

Problem	Possible Cause	Comments and suggestions
Poor DNA quality or yield	Deterioration or contamination of reagents.	Please ensure that the kit reagents are still in the effective using period before use. Discard any kit reagent that shows discoloration or evidence of microbial contamination.
	Kit stored under non-optimal conditions	Store kit at 15-25°C at all time after arrival. If either reagent or buffer precipitate upon shipping in cold weather or during long-term storage, dissolve precipitates by gently warming and stirring the solution. Please do not freeze the Reagent Cartridges.
	Insufficient sample input	DNA yield depends on the sample type and the number of nucleated cells in the sample. Please proportionally adjust the total input amount of sample to increase the DNA yield.
	Too much of elution buffer was used	The elution volume can be reduced proportionally.
	The eluate of final product(s) is not enough.	Please collect issue information and provide it to your Support Representative/Technical Support as soon as possible.
Clogging issue	Too much sample material was used.	Decrease the input amount of sample material or dilute your sample.
No results in downstream analysis	No signal/The PCR was inhibited.	Using appropriate controls for analysis. Check the positive control, negative control, water (NTC) and internal control to clarify the possible causes.
Instrument malfunction/abnormal sound	Abnormal consumables: 1. Deformed Filtered Tips 2. Deformed Reaction Chamber 3. Deformed Tip Holder	Please replace the batch with normal consumables.
	Abnormal action of instrument: 1. Inaccurate correction value 2. Spare parts or components damaged	Please collect issue information (videos and pictures) and provide it to your Support Representative/Technical Support as soon as possible to calibrate or replace any other damaged or worn parts.
*This table is helpful for solving common problem. If you need other technical support, please contact Zinexts team ( <a href="http://www.zinexts.com/index.php?lang=en">http://www.zinexts.com/index.php?lang=en</a> ) or your distributor.		

## Related Products








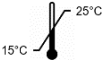



Product Name	Cat. no.
MagPurix <sup>®</sup> Blood DNA Extraction Kit 200	ZP02001
MagPurix <sup>®</sup> Blood DNA Extraction Kit 1200	ZP02002
MagPurix <sup>®</sup> Viral NA Extraction Kit	ZP02003
MagPurix <sup>®</sup> Tissue DNA Extraction Kit	ZP02004
MagPurix <sup>®</sup> Cultured Cell DNA Extraction Kit	ZP02005
MagPurix <sup>®</sup> Bacterial DNA Extraction Kit	ZP02006
MagPurix <sup>®</sup> HPV DNA Extraction Kit	ZP02007
MagPurix <sup>®</sup> TB DNA Extraction Kit	ZP02008
MagPurix <sup>®</sup> FFPE DNA Extraction Kit	ZP02009
MagPurix <sup>®</sup> Forensic DNA Extraction Kit	ZP02010
MagPurix <sup>®</sup> Viral Pathogen DNA Extraction Kit A	ZP02011
MagPurix <sup>®</sup> Viral Pathogen DNA Extraction Kit B	ZP02012
MagPurix <sup>®</sup> Viral RNA Extraction Kit	ZP02013
MagPurix <sup>®</sup> Plant DNA Extraction Kit	ZP02014
MagPurix <sup>®</sup> Total RNA Extraction Kit	ZP02015
MagPurix <sup>®</sup> Viral NA Extraction Kit LV	ZP02016
MagPurix <sup>®</sup> CFC DNA Extraction Kit	ZP02017
MagPurix <sup>®</sup> cfDNA Extraction Kit Plus	ZP02024
MagPurix <sup>®</sup> cfDNA Extraction Kit LV	ZP02025
MagPurix <sup>®</sup> Coronavirus RNA Extraction Kit	ZP02027
MagPurix <sup>®</sup> Urine cfDNA Extraction Kit	ZP02032
MagPurix <sup>®</sup> Plasma cfDNA RNA Extraction Kit	ZP02033

## References

- Tan SC *et al.* J Biomed Biotechnol. (2009)
- Collection, Transport, Preparation, and Storage of Specimens for Molecular Methods; Approved Guideline
- <http://www.cdc.gov/od/ohs/biosfty/biosfty.htm>
- <http://www.cdc.gov/od/ohs/pdffiles/DOTHazMat8-14-02.pdf>
- [http://www.who.int/csr/resources/publications/WHO\\_CDS\\_CSR\\_LYO\\_2004\\_9/en](http://www.who.int/csr/resources/publications/WHO_CDS_CSR_LYO_2004_9/en)
-

# Symbols

The following symbols are used on labels and in Instructions for Use (IFU), in compliance with EN ISO 15223-1 standard.

Symbol	Explanation
	CE mark
	For In Vitro Diagnostic Use
	Catalogue number
	Lot/Batch number
	Sufficient for [n] samples
	Instructions for Use
	Expiry date
	Storage temperature (15°C – 25°C)
	Manufacturer
	European Authorized Representative
	Caution

## Limited Product Warranty

Zinexts Life Science is committed to provide customers with high-quality products and services. Our goal is to ensure that every customer is 100% satisfied with our products and services. If you have any question or concerns, contact our Technical Support Representatives.

Zinexts Life Science guarantees the performance of all products according to the specifications stated in our product literature. The purchasers/users must determine the suitability of the product for their particular use. We reserve the right to change, alter, or modify any product to enhance its performance and design.

This warranty limits Zinexts Life Science Corporation's liability only to the cost of the product. No warranty is granted for products beyond their listed expiration date. No warranty is applicable unless all product components are stored and used in accordance with instructions.

## Revision History

Version	Date	Description
2.0	11. Aug. 2022	List of IVD symbols added
2.1	13.Oct. 2022	1. binding buffer change to 8 2. correct the centrifuge speed of leukocyte concentration